



TheraFit Premium Soft Goods

Shoulder Harness - Vest - Chest Strap - 4 Point Hip Belt - **Hip Belt** - Rear Pull Hip Belt - Center Pull Hip Belt



TheraFit Hip Belts

- Available metal buckle w/push button release, or plastic side release buckle
- Pads are attached to the straps
- Metal attaching clips, hinged collars, and cam lock buckles attach belt to chair

1" Straps, 4.25"L x 1.75"H Pads	Metal Buckle	Plastic Buckle
w/Metal Clips	.61001	65001
w/Cam Locks	.61002	65002
w/Metal Clips & Cam Locks	.61003	65003
w/Metal Clips & Hinged Collars	.61004	65004
w/Cam Locks & Hinged Collars	.61005	65005
w/Clips, Cams & Hinged Collars	.61006	65006

1" Straps, 5.25"L x 1.75"H Pads	Metal Buckle	Plastic Bkl
w/Metal Clips	.61101	65101
w/Cam Locks	.61102	65102
w/Metal Clips & Cam Locks	.61103	65103
w/Metal Clips & Hinged Collars	.61104	65104
w/Cam Locks & Hinged Collars	.61105	65105
w/Clips, Cams & Hinged Collars	.61106	65106

1" Straps, 6.25"L x 1.75"H Pads	Metal Buckle	Plastic Bkl
w/Metal Clips	.61201	65201
w/Cam Locks	.61202	65202
w/Metal Clips & Cam Locks	.61203	65203
w/Metal Clips & Hinged Collars	.61204	65204
w/Cam Locks & Hinged Collars	.61205	65205
w/Clips, Cams & Hinged Collars	.61206	65206

1.5" Straps, 5.25"L x 2.25"H Pads	Metal Buckle	Plastic Buckle
w/Metal Clips	.62001	65301
w/Cam Locks	.62002	65302
w/Metal Clips & Cam Locks	.62003	65303
w/Metal Clips & Hinged Collars	.62004	65304
w/Cam Locks & Hinged Collars	.62005	65305
w/Clips, Cams & Hinged Collars	.62006	65306

1.5" Straps, 7.25"L x 2.25"H Pads	Metal Buckle	Plastic Bkl
w/Metal Clips	.62101	65401
w/Cam Locks	.62102	65402
w/Metal Clips & Cam Locks	.62103	65403
w/Metal Clips & Hinged Collars	.62104	65404
w/Cam Locks & Hinged Collars	.62105	65405
w/Clips, Cams & Hinged Collars	.62106	65406

1.5" Straps, 9.25"L x 2.5"H Pads	Metal Buckle	Plastic Bkl
w/Metal Clips	.62201	65501
w/Cam Locks	.62202	65502
w/Metal Clips & Cam Locks	.62203	65503
w/Metal Clips & Hinged Collars	.62204	65504
w/Cam Locks & Hinged Collars	.62205	65505
w/Clips, Cams & Hinged Collars	.62206	65506

2" Straps, 9.25"L x 3"H Pads	Metal Buckle	Plastic Buckle
w/Metal Clips	.64001	65601
w/Cam Locks	.64002	65602
w/Metal Clips & Cam Locks	.64003	65603
w/Metal Clips & Hinged Collars	.64004	65604
w/Cam Locks & Hinged Collars	.64005	65605
w/Clips, Cams & Hinged Collars	.64006	65606

2" Straps, 11.25"L x 3"H Pads	Metal Buckle	Plastic Bkl
w/Metal Clips	.64101	65701
w/Cam Locks	.64102	65702
w/Metal Clips & Cam Locks	.64103	65703
w/Metal Clips & Hinged Collars	.64104	65704
w/Cam Locks & Hinged Collars	.64105	65705
w/Clips, Cams & Hinged Collars	.64106	65706



1" STRAPS

1.5" STRAPS

2" STRAPS

TheraFit Hip Belt Length	Lengths are taken from
1" Straps51"	end to end, straps fully
1 1/2" Straps58"	extended, without any of
2" Straps67"	the attaching hardware.

TheraFit Hip Belts w/Metal Buckle - Bulk Packs Of 20

- Choose a belt and choose the same size attaching hardware
- Give both part numbers when ordering

Belt Bulk Pack, 20 Pieces	Attaching Hardware Bulk Pack, 20 Sets
1" Straps, 4.25"L x 1.75"H Pads . . .61000	1" Metal Clips63101
1" Straps, 5.25"L x 1.75"H Pads . . .61100	1" Cam Locks63102
1" Straps, 6.25"L x 1.75"H Pads . . .61200	1" Metal Clips & Hinged Collars . .63104
	1" Cam Locks & Hinged Collars . .63105
1.5" Straps, 5.25"L x 2.25"H Pads .62000	1.5" Metal Clips63111
1.5" Straps, 7.25"L x 2.25"H Pads .62100	1.5" Cam Locks63112
1.5" Straps, 9.25"L x 2.5"H Pads . . .62200	1.5" Metal Clips & Hinged Collars .63114
	1.5" Cam Locks & Hinged Collars .63115
2" Straps, 9.25"L x 3"H Pads64000	2" Metal Clips63121
2" Straps, 11.25"L x 3"H Pads64100	2" Cam Locks63122
	2" Metal Clips & Hinged Collars . .63124
	2" Cam Locks & Hinged Collars . .63125



Attaching Hardware Information on page 50


HIP BELTS & CHEST STRAPS

*Parts included vary depending on which belt you have.
You may not have all parts shown here.*


METAL ATTACHING CLIPS

Depending on width of strap, thread strap as shown.

Metal Attaching Clip & 1" Or 1.5" Strap

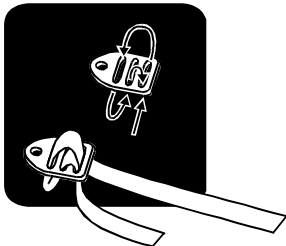


Metal Attaching Clip & 2" Strap

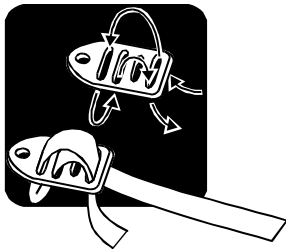


CAM LOCK BUCKLES

Thread strap through buckle and push top down to hold strap in place.

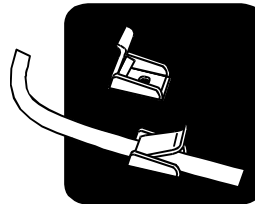


SLIDING DOUBLE-D CLIPS

Use Sliding Double-D Clips to loop straps around chair frame.

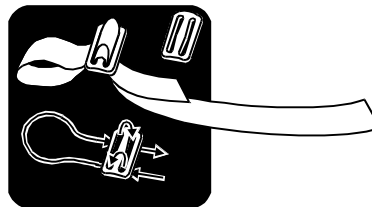


Use Sliding Double-D Clips with Metal Attaching Clips.

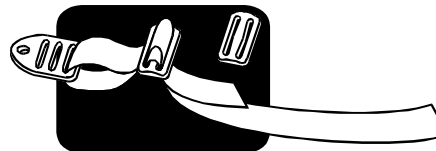

HINGED COLLARS

Attach Hinged Collars to chair frame. Attach Metal Attaching Clips or Cam Lock Buckles to Hinged Collars.

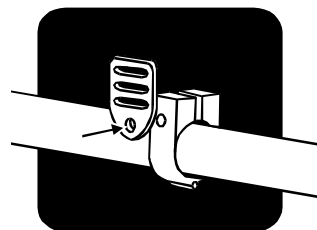

08756