


iontoPatch **EXTRA** STRENGTH®

Electronic Transdermal Iontophoretic Drug Delivery System

Stronger dose* for more aggressive yet safe treatment.

*120mA-min

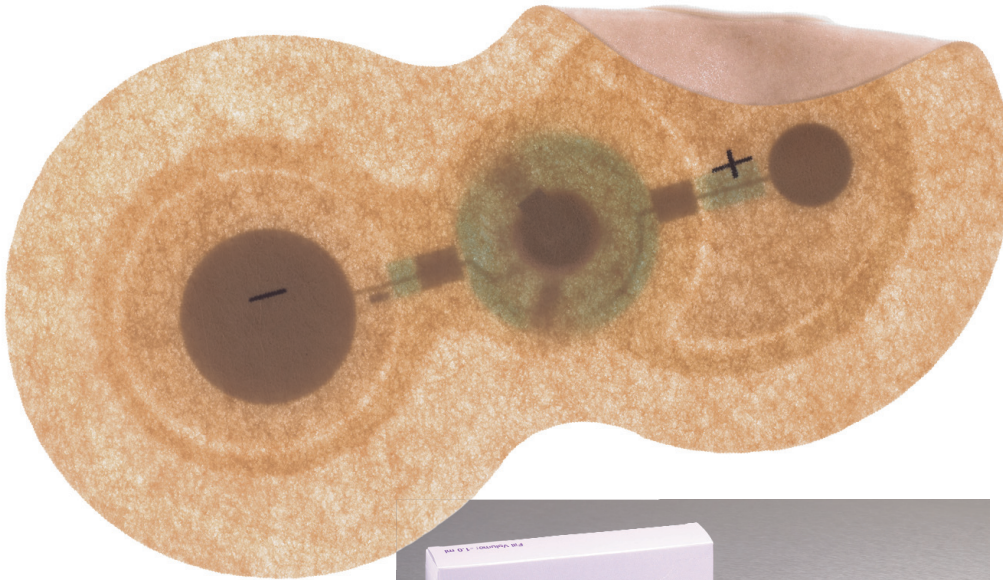

- **More patients improve faster to allow them to participate in full therapy treatment sooner**
- **Higher dose delivers more therapy in fewer treatments**
- **Safer for patients susceptible to skin irritation**
- **8 hours average patient wear time****
- **Treats most areas, including elbows, knees, wrists, feet and shoulders.**
- **Wrappable, compressible and braceable**
- **Extended time-released iontophoresis**

**Delivery time is dependent upon individual patient skin resistance.

iontoPatch Extra Strength® Specifications:

Active Area	1.9 in ² (12 cm ²)
Approximate Fill Volume	1.0 ml
Approximate Operating Current	0.25 mA
Dosage	120 mA-min
Average Wear-time	8 hours (Delivery time is dependent upon individual patient skin resistance.)
Skin Interface	Synthetic, non-woven absorbent material
Patch Adhesive	Breathable and highly flexible; latex-free
Conductive Element	Screen-printed electrode and lithium battery
Packaging	6 patches per pack, 6 saline vials
Product Number	D-0087-120

Caution: Federal Law (USA) restricts this device to use by or on the order of a physician.



Our complete iontoPatch family of products: iontoPatch® SP, iontoPatch® STAT®, iontoPatch® 80, iontoPatch® Extra Strength®

See the link to the iontoPatch Application Instruction Video located on the iontoPatch page of the Patterson Medical website at www.PattersonMedical.com

Patterson Medical product code:
 iontoPatch Extra Strength® 5509-36
 iontoPatch® 80 A939-65-0
 iontoPatch® STAT® 5640-91
 iontoPatch® SP 9226-40

PATTERSON
MEDICAL

TM-0101-0813 Travanti Medical® ©2013