

R3xm *recumbent cycle*

Cardio product solutions

BELIEVE IN BETTER HEALTH

Our R3xm Recumbent Cycle is specially designed to meet Med CE IIa certification standards, making it possible to incorporate low-impact cardio with a controlled range of motion early in your clients' continuum of care. A step-through design offers easy entry and exit, and almost everything about this cycle can be adjusted to fit your clients' unique needs. Pedal cranks are independently height-adjustable, making it possible to minimize the movement of the knee or hip with limited range of motion. Seat position can be adjusted with one hand, and the seat swivels from full right to full left to assist with ingress and egress. Both pedals and handles are thoughtfully designed for comfort and stability, and the easy-to-navigate console streamlines feedback to show you exactly what you need to see to help your clients down the road to wellness.




FEATURES	
Contact & Telemetric HR	Yes
Handlebar Design	Seat side / Front ergo bend
Integrated Arm Rests	No
Remote Buttons	Yes
Bottle Holder	Yes

CONSOLE	
Display Type	Dot-matrix LED with profile display
Display Feedback	Time of Day, Time Elapsed, Time Remaining, Total Program Time, Distance (Kilometers or Miles), Calories, Level, Speed, RPM, Heart Rate, METs, Watts, Workout Profile
Secondary Data Display	No
User-defined Multi-language Display	Accessible only through manager menu: English, German, French, Italian, Spanish, Dutch, Portuguese, Japanese, Swedish, Finnish, Turkish, Polish
Resistance Range	1-30
Workouts	Manual, Constant Watts, Fat Burn, Interval Training, Target Heart Rate, Fitness Test
CSAFE Ready-FitLinxx™ Certified	Yes
IPTV Compatible	No
Pro:Idiom Compatible	No
FitTouch™ Technology	No
Instant On Technology	Yes
One-button Start	Yes
On-the-fly Program Change	Yes
Integrated Vista Clear™ Television Technology	No
FITCONNEXION™ Ready	No
WiFi Enabled	Optional
Bluetooth Enabled	No
iPod® / iPhone® / iPad® Compatible	Charging only
Multimedia Playback	No
Personal Fan	No
Personal Trainer Portal Compatible	No
USB Port	Yes
Virtual Active™ Compatible	No
Asset Management Compatible	Yes
Workout Tracking Network Compatible	No
Web Connectivity	No
Facility Communication Portal/Calendar	No
RFID Compatible	No
Pause Function	Yes
Reading Rack	Yes



Dot-matrix with Profile Display

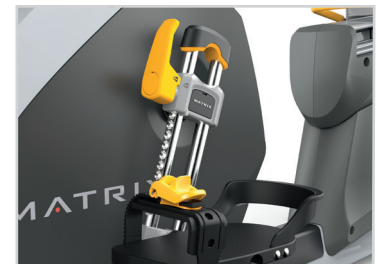

Height-adjustable Pedal Crank



Swivel Seat



One-hand Seat Adjustment

RESISTANCE SYSTEM	
Resistance System	Brushless generator

TECH SPECS	
Crank Design	User adjustable cranks
Power Requirements	Self-powered / powered 100v-240v – 50/60Hz AC
Minimum Watts	4 powered / 7 self-powered
Minimum RPM	25 powered / 25 self-powered
Assembled Dimensions	150.3 x 65.4 x 140 cm / 59.2" x 25.7" x 55.1"
Max User Weight	182 kg / 400 lbs.
Assembled Weight	87.5 kg / 192.9 lbs.
Shipping Weight	99.7 kg / 219 lbs.