

# The Authentic HIVAMAT™ Deep Oscillation



[hivamatus.com](http://hivamatus.com)



Personal



Evident Two-Channel

## HIVAMAT Deep Oscillation

by **PHYSIOMED**®  
*The Creators of Deep Oscillation*

**Hivamat Evident**  
In-Clinic Professional Use

**Hivamat Personal**  
In-Clinic & At-Home Use  
For Professionals and Patients

### **CROSS●●●LINK**

*Transfer protocols from Evident to Personal*




**Patented Technology**

*US Patent# 7,343,203*



**30+ years of research  
and clinical application**

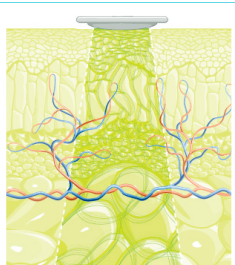

### **Non-Invasive & Deep-Acting Clinically Proven Effects:**

- Relieves Pain
- Reduces Secondary Lymphedema
- Relaxes Muscles
- Improves Lymphatic Circulation
- Promotes Mobility
- Activates Connective Tissue

### **Deep Oscillation**

Deep Oscillation is a unique, patented, proven, non-invasive, athermal therapy method that creates biological oscillations using electrostatic attraction and friction. This generates a gentle kneading effect deep within the tissue.

**Electrostatic Oscillations Permeate Up To 8 cm Deep**

## Product Features

### Frequency Range

Selection frequency from 5 - 250 Hz and 4 sweeping ranges

### One-Button Operation

Single data selector for fast on/off and simple to use interface

### Color Display

Multi-color LCD display  
**Evident:** Touch Screen

### Channels

**Personal:** One Channel  
**Evident:** Two Channels

### Power

**Personal:** 7.2v rechargeable lithium Ion  
**Evident:** 100 - 240 AC

### Dual Application

Apply treatment with applicator or vinyl gloves for device-assisted manual therapy

### Therapy Card

A single card activates treatment protocols and can be transferred from the Evident to Personal

**CROSS●●●LINK**

### Treatment Programs

Pre-Programmed parameters can be dynamically adjusted during treatment or create individual settings of frequency, frequency bands, duty cycle, duration of therapy and burst function


### Control Light

A blue light beneath the data selector indicates treatment in progress and intensity



### Anatomical Animations

Anatomical visual overview of treatment

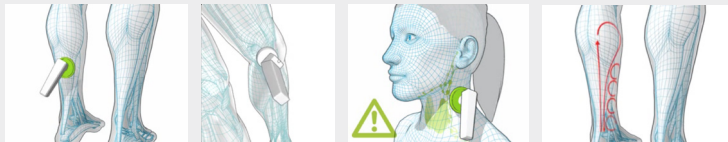

### Warranty

Hivamat Personal and Evident come standard with a 3-year warranty



#### Personal

In-Clinic & At-Home Use



#### Evident

In-Clinic Professional Use

### Special Evident Features



Swivel and tilt mount for optimum monitor alignment at all times



8" Color Monitor – Combination touch screen and dial control



100 Patient Database – Transfer programs to Personal unit

### Accessories Included

Applicator handhold (1)  
Oscillator head 5 cm (2)  
Oscillator head 9.5 cm (2)  
Hospital-Grade Mains Cable (Evident)  
Connection Cable Deep Oscillation (2)  
Electrode Connection Cable (1)

Patient Lead (Evident)  
Therapy card (1)  
Blank Therapy Cards (3 - Evident)  
Titanium neutral element (1)  
Box of vinyl gloves (1)  
Self-adhesive 3.2 cm snap electrode (1 pk. of 4)