

NAME: _____

PHONE: _____

SCHOOL/ORGANIZATION: _____

E-MAIL: _____

DEALER/REP: _____

SALES REP: _____ QUOTE: _____



STANDARD FEATURES:

- Load Capacity: For use with patients up to 600 lbs.
- 20.5"W x 30"D x 50-60"H (Assembled)
- White powder coat base
- Black upholstery
- Includes 1 Box (50) Sani-Cover® Fitted Disposable Face Rest Covers
- Epidural Positioning Device (EPD)
- 180° Adjustable face rest
- 6 Position adjustable arm rest
- Arm rest size: 16.5"W x 10"D
- Arm rest adjustment: 32-43" from the floor
- Interchangeable torso cushions (1 flat and 1 elevated)
- Torso cushion size: 8"W x 8"H
- Torso cushion adjustment: 5.5-10"
- EPD Stand (LDS)
- Positive locking clamp for EPD
- 7 Position adjustable foot rest
- Foot rest size: 20"W x 12"D
- Foot rest adjustment: 10.5-19.5" from the floor
- Medical grade profile locking casters

EPD Package 2

EPIDURAL POSITIONING DEVICE AND STAND



*EPD Package 2
(Epidural Positioning
Device + Stand)
with flat torso cushion*



*Elevated torso
cushion*



FRC 30-50



CE Patents:
US 6,758,447, B2
and D455,495

Rethink patient safety with the patented and original Epidural Positioning Device. This revolutionary device features unique torso positioning, exclusive to our design, that moves the back into an arched position for proper cervical alignment to position patients safely and comfortably during epidural procedures. Improper patient positioning during epidural anesthesia is the number-one cause of complications in the operating room.

The EPD Package 2 includes both the Epidural Positioning Device (EPD) and Labor Delivery Stand (LDS) for maximum versatility.

Epidural Positioning Device (EPD)

Patented in 2001, our Epidural Positioning Device is a landmark device designed by an anesthesiologist that positions patients correctly and comfortably by encouraging cervical, thoracic and lumbar flexion while maintaining a solid and stable position. Manually controlled chest, arm and head support ensures that patients of every body type and size can be accommodated and the back remains immobile during epidural placement. The EPD is suitable for use in operating rooms, labor and delivery suites, pain centers or anywhere epidural and spinal procedures are performed.

EPD Stand (LDS)


Designed to be used with the EPD, the EPD stand provides adjustable foot support to aid in positioning the patient and for patient comfort. It also provides an adjustable mount for the EPD and features a built-in clamp allowing for the same range of adjustment as an OR table.

EPD Package 2

EPIDURAL POSITIONING DEVICE AND STAND

STANDARD PACKAGE:

- EPD Package 2**
Epidural Positioning Device & EPD Stand
Qty _____

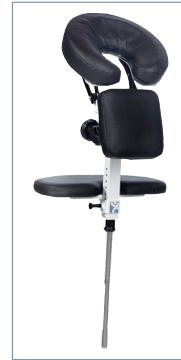

EPD - Epidural Positioning Device



LDS - EPD Stand

OPTIONS:

- FRC 30-50** Sani-Cover® Fitted Disposable Face Rest Covers (50 Covers per Box)
Qty _____
- FRC 30-50-CASE** Sani-Cover® Fitted Disposable Face Rest Covers (50 Covers per Box) Case of 18 Boxes (900)
Qty _____



Sani-Covers™ can be used on the face rest and the arm rest of the EPD