



Accommodator™ semi-rigid PRO

Comfort of a soft insole with the control of a semi-rigid orthotic


ACCOMMODATOR™ SEMI-RIGID PRO

Combine the comfort of a soft insole with the control of a semi-rigid orthotic. By stabilizing the foot, supporting the arch, and reducing pronation, these prefabricated orthotics offer the wearer immediate relief from heel and arch pain—right out of the box!

The base of this multi-layered insole is resilient, expanded EVA foam offering cushioning and shock absorption from heel to toe. Between the base and top cover is a structured, semi-rigid shell that cradles the heel and supports the arch, offering moderate pronation control. Unlike other semi-rigid orthotics, this shell has a Poron™-filled heel window strategically positioned beneath the plantar fascia origin, maximizing heel-strike protection. The entire insole is cushioned by a layer of breathable foam that wicks away moisture and limits perspiration. The foam is laminated to an antimicrobial-treated fabric that helps control odor and friction.

Available in Full-Length version or ¾-Length to fit a wider variety of shoes. Full-length insoles can be scissor-trimmed for a custom fit. Insoles can be easily modified through the utilization of posting strips.

- #0815-80-158 Full Size AA
- #0815-80-166 Full Size A
- #0815-80-174 Full Size B
- #0815-80-182 Full Size C
- #0815-80-190 Full Size D
- #0815-80-208 Full Size E
- #0815-80-216 Full Size F
- #0815-80-224 Full Size G
- #0815-80-232 Full Size H
- #0815-80-240 Full Size I/J
- #0815-80-257 Full Size K
- #0815-80-265 ¾ Size AA
- #0815-80-273 ¾ Size A
- #0815-80-281 ¾ Size B
- #0815-80-299 ¾ Size C
- #0815-80-307 ¾ Size D
- #0815-80-315 ¾ Size E
- #0815-80-349 ¾ Size F

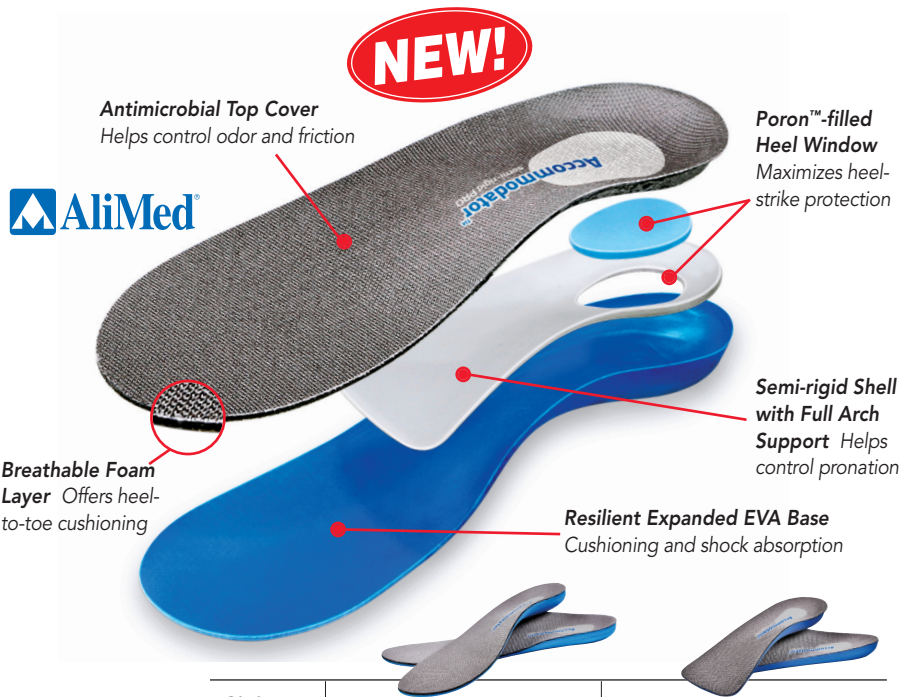

Sizing Insole Size	Full-Length		¾-Length	
	Wmns	Mens	Wmns	Mens
AA	5-5½	3-3½	6-6½	4-4½
A	6-6½	4-4½	7-8½	5-6½
B	7-7½	5-5½	9-10½	7-8½
C	8-8½	6-6½	11-11½	9-9½
D	9-9½	7-7½	12-12½	10-10½
E	10-10½	8-8½	13	11-11½
F	11-11½	9-9½	14-15	12-13½
G	12-12½	10-10½	—	—
H	13	11-11½	—	—
I/J	14-15	12-13	—	—
K	—	14	—	—

KEY BENEFITS

- ▶ Natural orthotic with full arch support and deep heel cup help control pronation
- ▶ Semi-rigid insert with poron-filled heel window absorb shock forces during heel strike
- ▶ Foam layer with antimicrobial-treated top cover minimizes friction and controls odor



www.pattersonmedical.com

Call toll free:
1-800-323-5547
 In Canada call 1-800-665-9200

Fax toll free:
1-800-547-4333
 In Canada call 1-877-566-8611