

## AliMed® BreakAway Lap Cushion

#76808, #70065, #71694, #71 0092 16" x 18" • #71349 24" x 26"


### Directions for Applications, Use & Care

#### Application

The BreakAway Lap Cushion is a convenient, non-restraint device used to safely secure the wheelchair occupant in the seated position. The cushion is designed to be secured in place with clip on one side of the cushion, and releasable hook and loop fastener on the other side. When the patient wants to stand, the hook and loop fastener will release. The patient does not need to manually remove the pad. The Lap Cushion, with a water-resistant surface, can also be used as a work surface.

#### Installation of Lap Cushion

1. Lock the wheelchair brakes and rotate the footrests to the side of the wheelchair.
2. Seat the occupant as far back in the wheelchair as possible.
3. With wheelchair occupant seated, position the lap cushion between the top of their legs and the underside of the wheelchair armrests. Be sure to face the cutout side of the cushion toward the occupant. (see Fig. 1)
4. Wrap the webbing with the clip around the front vertical tube of the wheelchair armrest and clip to D-ring.
5. Secure the hook and loop fastener around the other vertical tube. (see Fig. 2)
6. Make sure the magnets make a good contact. (see Fig. 3)
7. Route sensor cable to Alarm unit. Connect sensor cable to Alarm unit.
8. Test

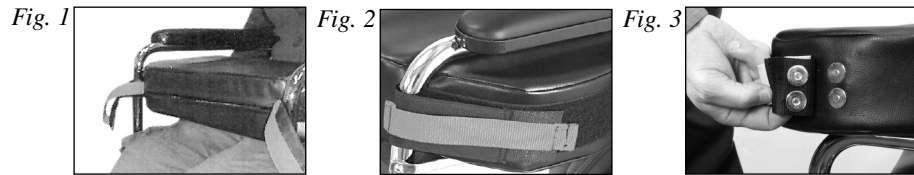

Care

Clean cushion with mild detergent and/or disinfectant. Use a damp cloth or sponge. Do not launder or submerge in any liquid.

#### To Reorder

#MJ76808 BreakAway Lap Cushion, w/o Alarm 16" x 18"W  
#MJ71349 BreakAway Lap Cushion w/o Alarm, 24" x 26"W

#### Other Alarm and Safety Products

#MJ75188 BreakAway Lap Cushion w/Basic Alarm  
#MJ77643 Floor Sensor Mat w/Alarm  
#MJ77655 2" Folding Fall Pad w/Alarm  
#MJ77656 4" Folding Fall Pad w/Alarm  
#MJ70065 Cushion w/TR2 Alarm  
#MJ2143 Wireless Bed Alarm System  
#MJ2144 Wireless Chair Alarm System  
#MJ77067 Bed Sensor System  
#MJ77062 Chair Sensor System  
#MJ75049 Sensor Cushion Alarm  
#MJ71694 Cushion w/ IQ Duo Alarm  
#MJ710092 w/Voice Alarm

Distributed by:

 **AliMed® inc.**  
Dedham, MA 02026-9135

www.AliMed.com  
#5938  
©2008 AliMed® inc. All rights reserved.

## AliMed® BreakAway Lap Cushion

#76808, #70065, #71694, #71 0092 16" x 18" • #71349 24" x 26"


### Directions for Applications, Use & Care

#### Application

The BreakAway Lap Cushion is a convenient, non-restraint device used to safely secure the wheelchair occupant in the seated position. The cushion is designed to be secured in place with clip on one side of the cushion, and releasable hook and loop fastener on the other side. When the patient wants to stand, the hook and loop fastener will release. The patient does not need to manually remove the pad. The Lap Cushion, with a water-resistant surface, can also be used as a work surface.

#### Installation of Lap Cushion

1. Lock the wheelchair brakes and rotate the footrests to the side of the wheelchair.
2. Seat the occupant as far back in the wheelchair as possible.
3. With wheelchair occupant seated, position the lap cushion between the top of their legs and the underside of the wheelchair armrests. Be sure to face the cutout side of the cushion toward the occupant. (see Fig. 1)
4. Wrap the webbing with the clip around the front vertical tube of the wheelchair armrest and clip to D-ring.
5. Secure the hook and loop fastener around the other vertical tube. (see Fig. 2)
6. Make sure the magnets make a good contact. (see Fig. 3)
7. Route sensor cable to Alarm unit. Connect sensor cable to Alarm unit.
8. Test

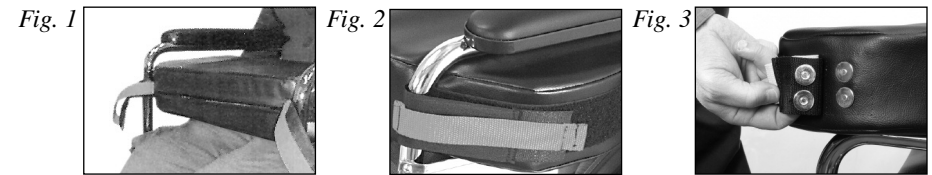

Care

Clean cushion with mild detergent and/or disinfectant. Use a damp cloth or sponge. Do not launder or submerge in any liquid.

#### To Reorder

#MJ76808 BreakAway Lap Cushion, w/o Alarm 16" x 18"W  
#MJ71349 BreakAway Lap Cushion w/o Alarm, 24" x 26"W

#### Other Alarm and Safety Products

#MJ75188 BreakAway Lap Cushion w/Basic Alarm  
#MJ77643 Floor Sensor Mat w/Alarm  
#MJ77655 2" Folding Fall Pad w/Alarm  
#MJ77656 4" Folding Fall Pad w/Alarm  
#MJ70065 Cushion w/TR2 Alarm  
#MJ2143 Wireless Bed Alarm System  
#MJ2144 Wireless Chair Alarm System  
#MJ77067 Bed Sensor System  
#MJ77062 Chair Sensor System  
#MJ75049 Sensor Cushion Alarm  
#MJ71694 Cushion w/ IQ Duo Alarm  
#MJ710092 w/Voice Alarm

Distributed by:

 **AliMed® inc.**  
Dedham, MA 02026-9135

www.AliMed.com  
#5938  
©2008 AliMed® inc. All rights reserved.